



### Key Features

- Stores completely flat
- Radio translucent
- CT and MRI compatible
- Sizing line for easy and accurate sizing
- Standard size markings on both Adult and Mini
- Chin flip feature to ease airway management
- Nasal cannula tubing holder support
- Facilitates intubation
- Latex free
- Disposable
- When shipped, the default position of the the Adult Ace collar is in the Neckless position, and the default position of the Mini Ace is the Infant position

## Ambu<sup>®</sup> Perfit ACE<sup>®</sup> & Mini Perfit ACE<sup>®</sup>

### 28 sizes in two collars

The Ambu Perfit ACE is a unique adjustable collar. Combined with the Mini Perfit ACE these two collars offer an incredible 28 sizes ensuring that patients can now have the most effective and accurate fit ever.

Made from radio translucent materials, the collar is packaged completely flat and has a unique chin-flip feature to ease airway management. Ambu Perfit ACE is a disposable adult adjustable collar for immobilizing the cervical column.

Ambu Perfit Ace is an adjustable collar designed to fit any adult with 16 different settings ranging from Neckless to Tall sizes. The Ambu Perfit ACE incorporates latches that allow the rescuer to size the collar exactly to the victim's neck size. Two safety lock buttons hold the collar in place when adjusted.

Ambu Mini Perfit ACE is an adjustable collar designed to fit children and small adults, with 12 different settings to provide an accurate fit.

### MATERIALS – AMBU PERFIT ACE

White/light grey plastic	Polypropylene
Foam	Polyethylene

### MATERIALS – AMBU MINI PERFIT ACE

Yellow plastic	Polypropylene
Foam	Polyethylene

IMMOBILIZATION

## SPECIFICATIONS – AMBU PERFIT ACE

Adjustable sizes	Neckless, Short, Regular and Tall
Number of settings	16 adjustments
Dimensions	554 x 215 x 14 mm (21.8" x 8.5" x 0.6")
Weight	200 g (7 ounces)
Color	White/light grey
Translucency	CT and MRI compatible



Ambu® Perfit ACE®

## SPECIFICATIONS – AMBU MINI PERFIT ACE

Adjustable sizes	Infant, Pediatric and Small Adult
Number of settings	12 adjustments
Dimensions	447 x 188 x 14 mm (17.6" x 7.4" x 0.6")
Weight	100 g (3.5 ounces)
Color	Yellow/light grey
Translucency	CT and MRI compatible




Ambu® Mini Perfit ACE®

## ORDERING INFORMATION

	CAT. NO.
Perfit ACE Extrication Collar, 30/case	000 281 000
Perfit ACE packaged with Head Wedge CID, 30/Case	000 281 151
Mini Perfit ACE Extrication Collar, 30/Case	000 281 106
Mini Perfit ACE packaged with Head Wedge CID, 30/Case	000 281 161

Available packaged with the Ambu® Head Wedge®.

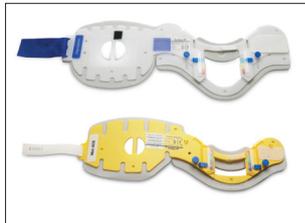

Ambu® Perfit ACE®  
Ambu® Mini Perfit ACE®



Ambu® Head Wedge®



### Ambu A/S

Baltorpbakken 13  
DK-2750 Ballerup  
Denmark  
T +45 72 25 20 00  
F +45 72 25 20 50  
www.ambu.com

IMMOBILIZATION

