

# SINGLE-USE LIKE YOU'VE NEVER SEEN IT BEFORE

THE ULTIMATE TOOLBOX FOR BRONCHOSCOPY




Ambu® aScope™ 5 Broncho

**Ambu**

# ONE SOLUTION

# TO EQUIP YOU FOR A BROAD RANGE OF PROCEDURES

## aScope 5 Broncho

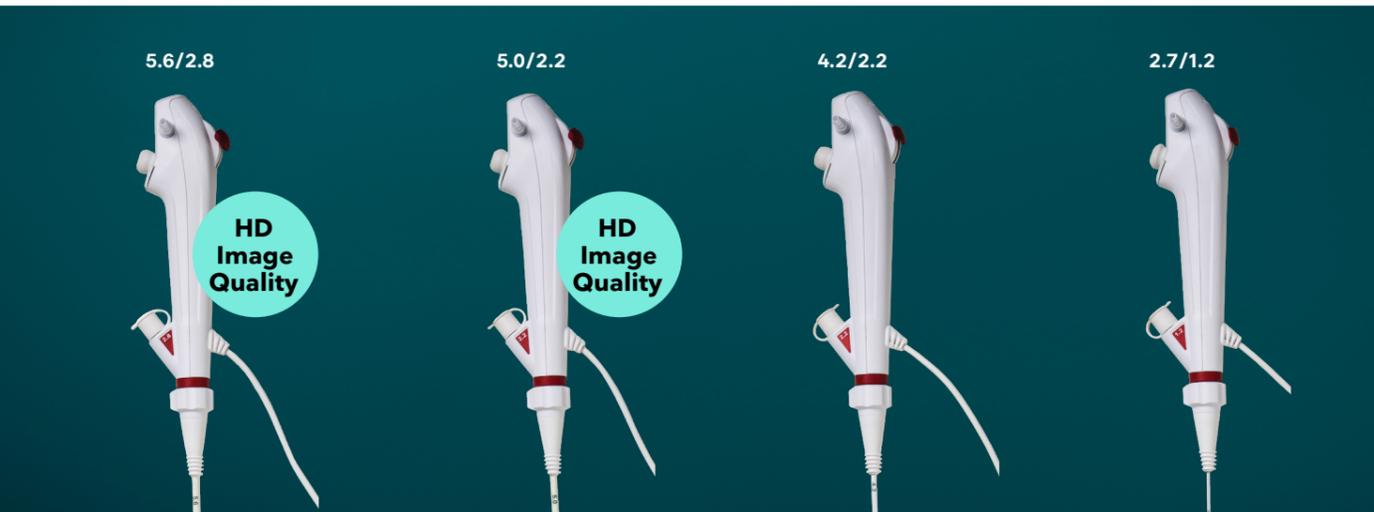

## aScope 5 Broncho Sampler Sets



### DIAGNOSTIC & THERAPEUTIC PROCEDURES

The entire aScope 5 Broncho range offers excellent manoeuvrability with bending angles of 195°/195° on the 5.6/2.8 and 5.0/2.2 sizes and 210°/210° on the 4.2/2.2 and 2.7/1.2 sizes; and even when working with tools, they continue to offer high degrees of flexion.

aScope 5 Broncho is indicated for use with a broad spectrum of endotherapy instruments, including high frequency electrosurgical tools\*. And the 4.2/2.2 size in particular offers a unique opportunity to take advantage of a regular-sized working channel in a thin bronchoscope.

All of this, combined with their bright and clear image, make them suitable for a wide range of diagnostic and therapeutic procedures.

### LARGER & SMALLER AIRWAYS


Whilst the two larger sizes of aScope 5 Broncho are ideal for the central airways, the thin and ultra-thin sizes are designed with smaller outer-diameters for navigating the outer peripheries of the lungs.

The two smaller sizes may also be suitable for those with smaller anatomies, with the 2.7/1.2 in particular being compatible with ET tubes as small as 3.5mm, whilst still offering an accessible working channel. Or they can simply be used to make a procedure more comfortable for patients.

4.2/2.2 & 2.7/1.2



5.6/2.8 & 5.0/2.2

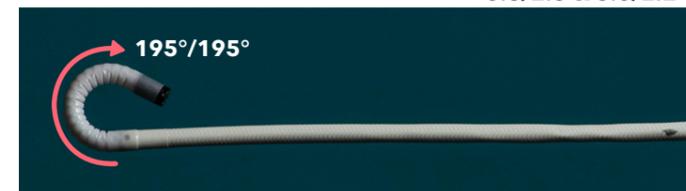

### BALS & BWS

The single-use aScope 5 Broncho Sampler Set is a unique, integrated sampling solution that protects the sample, simplifies workflow and improves safety.

It includes two 60 ml sample containers to address the need for larger samples in BAL and BW procedures.



#### Protects the sample from start to finish

Closed-loop, secure system reduces the risk of sample contamination and sample loss common with Lukens traps, to maximise quality and yield.

#### Simple to operate without assistance

Integrates seamlessly with the aScope 5 Broncho endoscope.

Allows switching between suction and sampling at the turn of a knob, and without the need to switch tubing.

Operators can easily change and close containers without the need for support staff.

#### Sterile and ready to go

Single-use solution is sterile straight from the package.

Always available when you need it.

Clicks right onto the aScope 5 Broncho endoscope.

\*aScope 5 Broncho 2.7/1.2 and 4.2/2.2 are not compatible with high-frequency tools

# PREPARING YOU FOR ALL EVENTUALITIES

## 100% UPTIME

Our single-use aScope 5 Broncho bronchoscopes are always available off the shelf. They require no reprocessing, repairs or servicing, meaning you'll experience no delays to your procedures as a result of scopes not being ready.

## REACT TO THE NEEDS OF YOUR CASES

We know your procedures don't always go to plan, and you may need a different size bronchoscope than you'd originally anticipated. With aScope 5 Broncho, you can quickly and easily react to the needs of your cases. With four variants to choose from and all being ready straight from their packs, you'll have the flexibility to switch to meet the evolving needs of your procedure without causing additional reprocessing burden or impacting scope availability for subsequent cases.

## FILL GAPS IN YOUR FLEET

Do you have all the bronchoscope sizes you need to tackle all the procedures you may be presented with? Are there are gaps in your bronchoscope fleet due to budgetary constraints, either because scopes haven't been purchased yet or they're broken and haven't been replaced? With aScope 5 Broncho, you can fill those gaps and maximise your capabilities, by choosing from four sizes, all available with minimal investment costs.



### STERILITY AND EFFICIENCY OF THE SINGLE-USE CONCEPT

- **Sterile:** No risk of patient-to-patient contamination
- **Always available:** No reprocessing or repair bottlenecks
- **Consistent performance:** No wear and tear affecting the bending performance or imaging



### THE LEVEL OF QUALITY NEEDED IN THE ENDOSCOPY UNIT

- **High performance:** Excellent manoeuvrability and high-quality imaging
- **Greater flexibility:** A variety of readily available bronchoscopes in 4 different sizes
- **A unique sampling solution:** The aScope 5 Broncho Sampler Set simplifies BAL procedures (see previous page)

# INTEGRATION

## WITH OUR LATEST SYSTEMS

### aBox 2 FULL-HD DISPLAYING & PROCESSING UNIT

aScope 5 Broncho bronchoscopes are compatible with the full-HD Ambu® aBox™ 2 displaying and processing unit. Its flexible and innovative design provides enhanced processing to deliver the outstanding image quality that's needed for the complex and demanding workflow of the endoscopy unit. Its compact design comprises of a light source, recorder and navigation software, and it can be connected to external monitors for a large-screen view.

### aCart Compact TROLLEY

The aCart workstations are space saving and provide convenient storage for Ambu displaying and processing units and endoscopes.



# LEADING THE WAY

## WITH SUSTAINABLE PRACTICES & TARGETS

From raw material extractions to disposal, the choices we make throughout the life cycle of our products not only affect our own operations, but also our suppliers and customers.

Our plastic and product strategy focuses on the areas where we believe we can make a difference in product design, production and disposal after end-of-use.

	SHORT TERM	MEDIUM TERM	LONG TERM
DESIGN	<b>Identify areas for improvement</b> Mature our capabilities where possible	<b>Reduce material use</b> <b>Eliminate use of PCV</b>	<b>Use sustainable materials</b> <b>Design for recycling</b>
PRODUCTION	<b>Operation Clean Sweep® program</b>	<b>Reduce unnecessary items in product packs</b>	<b>Waste to be reused or recycled</b>
END-OF-USE	<b>Plastic Bank® recycling partnership</b> <b>Take-back program with incineration &amp; energy recovery (US &amp; Italy)</b>	<b>Take-back program with recycling (other countries)</b>	<b>Recycling at scale</b>



Visit [www.ambu.com/sustainability](http://www.ambu.com/sustainability) or scan the QR code to read more about our targets.

Ambu® aScope™ 5 Broncho

# SINGLE-USE LIKE YOU'VE NEVER SEEN IT BEFORE

A new era for single-use bronchoscopy

## aScope™ 5 Broncho flexible single-use bronchoscope

Item code	Product	Quantity
622001000	Ambu® aScope™ 5 Broncho HD 5.6/2.8	5
621001000	Ambu® aScope™ 5 Broncho HD 5.0/2.2	5
620001000	Ambu® aScope™ 5 Broncho 4.2/2.2	5
624001000	Ambu® aScope™ 5 Broncho 2.7/1.2	5

## aScope™ 5 Broncho Sampler Set flexible single-use bronchoscope and specimen sampling system

Item code	Product	Quantity
622002000	Ambu® aScope™ 5 Broncho HD 5.6/2.8 Sampler Set	5
621002000	Ambu® aScope™ 5 Broncho HD 5.0/2.2 Sampler Set	5
620002000	Ambu® aScope™ 5 Broncho 4.2/2.2 Sampler Set	5

## Single-use sample container

Item code	Product	Quantity
480 500 000	Ambu® aScope™ BronchoSampler 60 SC	10

## Displaying and processing units

Item code	Product	Quantity
505001000EUC2	Ambu® aBox™ 2	1

## aCart™

Item code	Product	Quantity
9.713.800	Ambu® aCart™ Compact KU	1
3.000.991	Ambu® aCart™ Compact Vario Mount HA	1

**Ambu** FOREVER FORWARD

Ambu Limited  
Incubator 2  
Alconbury Weald Enterprise Campus  
Alconbury Weald  
Cambridgeshire  
PE28 4XA  
T: 01480 498403  
ambu.co.uk